


How to Enter results for non-DOCKED devices (POCT Manual Result) for Nurses

This tip sheet explains how to manually enter Point-of-Care Testing (POCT) results for non-docked devices in Epic using the Enter/Edit Results activity, highlighting that these results remain within Epic and do not interface to Cerner. This functionality is applicable to the following orders: Below are several Point-of-Care Testing (POCT) Manual Result tests available in Epic. Please note that this list is subject to change over time as new tests are created or existing ones are inactivated. When searching for a POCT test, use the search term 'POCT manual'

POCT Urinalysis Manual Dipstick Manual Result	POCT iSTAT ACTK Abbott Manual Result
POCT glucose Accucheck Manual Result	POCT iSTAT Creatinine Abbott Manual Result
POCT HIV-1 AND HIV-2 Oraquick HIV Manual Result	POCT nitrazine paper eye test Manual Result
POCT Urinalysis Clinitek Dipstick 10 Manual Result	POCT nitrazine paper test Manual Result
POCT Amniotest Manual Result	POCT Pregnancy Quick Vue Manual Result
POCT Clinitek Urine Pregnancy Manual Result	POCT PT/INR i-STAT PT Plus Manual Result
POCT fern test, vaginal fluid Manual Result	POCT Qualitative Urine Pregnancy Manual Result
POCT HDL/GLU/CHOL Cholestech Manual Result	POCT TC/HDL/TRIG Cholestech Manual Result
POCT Influenza A/B Quick Vue Manual Result	POCT Urinalysis Clinitek Dipstick 2 Manual Result

Enter/Edit Results

1. From the More Activities Drop-Down Icon, click the Enter/Edit Results Activity.

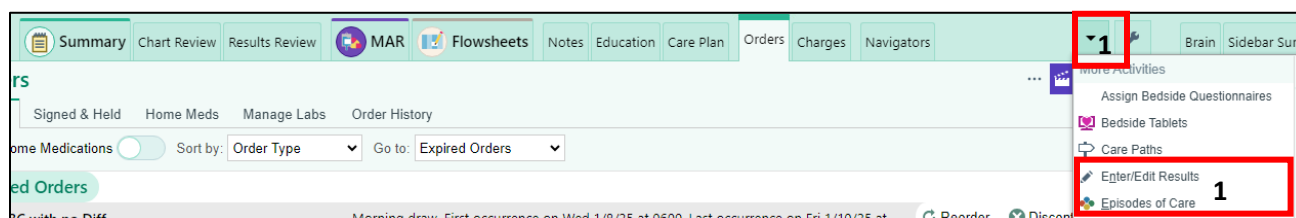

- From the Enter/Edit Results Activity, click on the necessary POCT manual test to result.
2. Enter specimen information.

Enter/Edit Results

☒ Preview
 ☒ Refresh
 ☒ Ext Result
 ☒ Edit
 ☒ Cancel Order
 ☒ Reflex Order

☒ Filters
 ☐ Outstanding Orders
 ☐ Default Types

Date/Time	Test	Auth. MD	Order Status	Result Status
Today				
Today at 18:24	POCT Glucose Accucheck Manual Result	Angelo J Acquista, MD	Active	

Searched through 7/9/2025 [Search Further](#)

☒ Scan
 ☒ Media Manager

POCT Glucose Accucheck Manual Result

Specimen

Type: Blood

Collected by: 3

Collection date: 3

Collection time:

☐ No collection information available

Resulting Lab

Lab name: LENOX HILL HOSPITAL 3

Technician: 3

Providers

Billing: ACQUISTA, ANGELO J

Resulting:

1. Enter specimen information. Collected by, Collection Date, Collection time and Technician are required fields. Resulting Lab, Lab Name will default to the facility the Patient is located in. NOTE: If the patient is in an outpatient location, the Lab Name will have to be manually entered.
2. Enter all the required result information for the test.

POCT Gastric Hemaprompt Manual Result

☐ No collection information available

Component	Value	Flags	Low	High	Ref Range	Units	Comment
POCT FOB CARD LOT#	2232						
POCT FOB CARD EXP DATE	8/7/2025				mm/dd/yy...		
POCT FOB DEVELOPER LOT#	3233						
POCT FOB DEVELOPER EXP DATE	8/7/2025				mm/dd/yy...		
POCT INTERNAL QC	Abnormal	Abnormal			Please i...		
POCT FOB CARD RESULT 1					Negative		
POCT FOB CARD RESULT 2					Negative		
POCT FOB CARD RESULT 3					Negative		

Result

Result date: 8/7/2025

Result time: 10:12 AM

Abnormal: Abnormal

Status: Final

- a. Lot # can be manually entered or scanned.
- b. If value of Fail is selected in Internal QC, the system will default Abnormal in the Flags column as well as the Abnormal field under Results. **NOTE:** When correcting Results, you will need to manually clear the Abnormal flags, they will not clear automatically.

- Result Status defaults to final for all Enter/Edit Results. Select Accept on the bottom right.

Enter/Edit Results from the Brain

The Inpatient RN can also perform Enter/Edit Results from the Brain.

- Click on the Go to Enter/Edit results in the Time Column
- Click on Enter/Edit Results

The screenshot shows the 'Brain' interface with patient information for 'One, Lijvpatient'. A task is listed at 1605: 'Go to Enter/Edit Results for POCT Glucose Accucheck Manual Result'. A red box labeled '1' highlights this task. Another red box labeled '2' highlights the 'Enter/Edit Results' link in the task details.

- The Enter/ Edit Result screen appears. The same steps will apply that are listed above for steps 3, 4, 4b and 5.

The 'Enter/Edit Results' screen displays the following information:

- Patient:** One, Lijvpatient, 14 y.o. / M, CSN: 11000007793, Code: Assume Full.
- Specimen:** Type: Blood, Collected by: [Redacted], Collection date: [Redacted], Collection time: [Redacted].
- Resulting Lab:** Lab name: LONG ISLAND JEWISH VALLEY STREAM, Technician: [Redacted].
- Providers:** Billing: CASE, JOSHUA DANIEL, Resulting: [Redacted].
- Components Table:**

Component	Value	Flags	Low	High	Ref Range	Units	Comment
POCT STRIP LOT # [3048800275]		!					
POCT STRIP LOT EXP DATE [30488...]		!			mm/dd/y...		
POCT EXTERNAL ACCEPTABLE QC...					Please i...		
POCT GLUCOSE ACCUCHECK [304...]		!	70	99		mg/dL	

Result: Result date: 8/8/2025, Result time: 10:32 AM, Abnormal: [Redacted], Status: Final.

Another way to Enter/Edit result

1. Click on the Labs icon from the brain

One, Lijvspatient
14 y.o. / M
CSN: 110000007793
Code: Assume Full

163-D

113

Orders

Labs 1

PRN

A S D

Attending Physician Inpatie...
Primary Problem
Burn (Principal Hospital Pr...
Allergies
Honey Bee Venom, Ging...

Infection
Bed bugs
Isolation
Cont
Notifications

2. Click on Enter Results

Lab Overview for One, Lijvspatient, 14 y.o. male, DOB: 3/4/2011, MRN: 10101465, Room: LIJVS 163, Bed: 163-D

Patient not scanned | Blood Collection: Lab

Available to Collect

Look Back: 24 Hours Show All Show: My Labs (1) All Labs (1)

Due	Priority	Order Name	Collect	Type	Source	Container	Actions
Yesterday 16:05	Routine	POCT Glucose Accucheck Manual Result	—	Blood	Blood, Capill...	—	Enter Results 2

3. The Enter/ Edit Result screen appears. The same steps will apply that are listed above for steps 3, 4, 4b and 5.

POCT Glucose Accucheck Manual Result

Specimen

Type: Blood

Collected by: [Redacted]

Collection date: [Redacted]

Collection time: [Redacted]

☐ No collection information available

Resulting Lab

Lab name: LONG ISLAND JEWISH VALLEY STREAM

Technician: [Redacted]

Providers

Billing: CASE, JOSHUA DANIEL

Resulting: [Redacted]

Cc Messages

+ PCP 1 + Care Team 2 + Treatment Team 3

+ Free Text 4 + My List 5 + Remove All 0

Add recipients

Send Results: No Yes Cc list only

Result

Result date: 8/8/2025

Result time: 11:15 AM

Abnormal: [Redacted]

Status: Final

Component	Value	Flags	Low	High	Ref Range	Units	Comment
POCT STRIP LOT # [3048800275]	[Redacted]	[Redacted]					
POCT STRIP LOT EXP DATE [30488...]	[Redacted]	[Redacted]			mm/dd/y...		
POCT EXTERNAL ACCEPTABLE QC...	[Redacted]	[Redacted]			Please i...		
POCT GLUCOSE ACCUCHECK [304...]	[Redacted]	[Redacted]	70	99		mg/dL	

