


How to Enter results for non-DOCKED devices (POCT Manual Result)

This tip sheet explains how to manually enter Point-of-Care Testing (POCT) results for non-docked devices in Epic using the Enter/Edit Results activity, highlighting that these results remain within Epic and do not interface to Cerner. This functionality is applicable to the following orders: Below are several Point-of-Care Testing (POCT) Manual Result tests available in Epic. Please note that this list is subject to change over time as new tests are created or existing ones are inactivated. When searching for a POCT test, use the search term 'POCT manual'

POCT Urinalysis Manual Dipstick Manual Result	POCT iSTAT ACTK Abbott Manual Result
POCT glucose Accucheck Manual Result	POCT iSTAT Creatinine Abbott Manual Result
POCT HIV-1 AND HIV-2 Oraquick HIV Manual Result	POCT nitrazine paper eye test Manual Result
POCT Urinalysis Clinitek Dipstick 10 Manual Result	POCT nitrazine paper test Manual Result
POCT Amniotest Manual Result	POCT Pregnancy Quick Vue Manual Result
POCT Clinitek Urine Pregnancy Manual Result	POCT PT/INR i-STAT PT Plus Manual Result
POCT fern test, vaginal fluid Manual Result	POCT Qualitative Urine Pregnancy Manual Result
POCT HDL/GLU/CHOL Cholestec Manual Result	POCT TC/HDL/TRIG Cholestec Manual Result
POCT Influenza A/B Quick Vue Manual Result	POCT Urinalysis Clinitek Dipstick 2 Manual Result

Enter/Edit Results

1. From the More Activities Drop-Down Icon, click the Enter/Edit Results Activity.

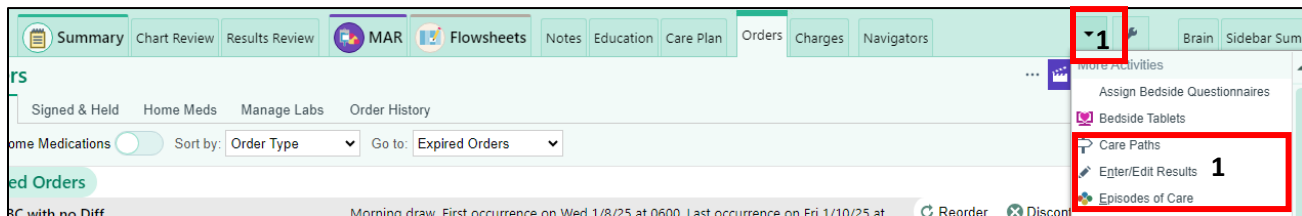

2. From the Enter/Edit Results Activity, click on the necessary POCT manual test to result.

Enter/Edit Results

☒ Preview ☒ Refresh ☒ Ext Result ☒ Edit ☒ Cancel Order ☒ Reflex Order

☒ Filters ☐ Outstanding Orders ☐ Default Types

Date/Time	Test	Auth. MD	Order Status	Result Status
Today				
Today at 18:24	POCT Glucose Accucheck Manual Result	Angelo J Acquista, MD	Active	

Searched through 7/9/2025 [Search Further](#)

☒ Scan ☒ Media Manager

POCT Glucose Accucheck Manual Result

Specimen		Resulting Lab	
Type: Blood		Lab name: LENOX HILL HOSPITAL	
Collected by:		Technician:	
Collection date:			
Collection time:			
☐ No collection information available		Billing: ACQUISTA, ANGELO J	
		Resulting:	

- Enter specimen information. Collected by, Collection Date, Collection time and Technician are required fields. Resulting Lab, Lab Name will default to the facility the Patient is located in. NOTE: If the patient is in an outpatient location, the Lab Name will have to be manually entered. If a provider is entering the Enter/Edit Results, they can enter their name in the Resulting Provider field as well as the Technician field.
- Enter all the required result information for the test.

POCT Gastric Hemaprompt Manual Result

☐ No collection information available		Resulting:																																									
<table border="1"> <thead> <tr> <th>Components</th><th>Sensitivities</th><th>Narrative</th><th>Impression</th></tr> </thead> <tbody> <tr> <td>Component</td><td>Value</td><td>Flags</td><td>Low High Ref Range Units Comment</td></tr> <tr> <td>POCT FOB CARD LOT#</td><td>2232</td><td></td><td></td></tr> <tr> <td>POCT FOB CARD EXP DATE</td><td>8/7/2025</td><td></td><td>mm/dd/y...</td></tr> <tr> <td>POCT FOB DEVELOPER LOT#</td><td>3233</td><td></td><td></td></tr> <tr> <td>POCT FOB DEVELOPER EXP DATE</td><td>8/7/2025</td><td></td><td>mm/dd/y...</td></tr> <tr> <td>POCT INTERNAL QC</td><td>Fail</td><td>Abnormal</td><td>Please i...</td></tr> <tr> <td>POCT FOB CARD RESULT 1</td><td></td><td>Negative</td><td></td></tr> <tr> <td>POCT FOB CARD RESULT 2</td><td></td><td>Negative</td><td></td></tr> <tr> <td>POCT FOB CARD RESULT 3</td><td></td><td>Negative</td><td></td></tr> </tbody> </table>		Components	Sensitivities	Narrative	Impression	Component	Value	Flags	Low High Ref Range Units Comment	POCT FOB CARD LOT#	2232			POCT FOB CARD EXP DATE	8/7/2025		mm/dd/y...	POCT FOB DEVELOPER LOT#	3233			POCT FOB DEVELOPER EXP DATE	8/7/2025		mm/dd/y...	POCT INTERNAL QC	Fail	Abnormal	Please i...	POCT FOB CARD RESULT 1		Negative		POCT FOB CARD RESULT 2		Negative		POCT FOB CARD RESULT 3		Negative		Result Result date: 8/7/2025 Result time: 10:12 AM Abnormal: Abnormal Status: Final	
Components	Sensitivities	Narrative	Impression																																								
Component	Value	Flags	Low High Ref Range Units Comment																																								
POCT FOB CARD LOT#	2232																																										
POCT FOB CARD EXP DATE	8/7/2025		mm/dd/y...																																								
POCT FOB DEVELOPER LOT#	3233																																										
POCT FOB DEVELOPER EXP DATE	8/7/2025		mm/dd/y...																																								
POCT INTERNAL QC	Fail	Abnormal	Please i...																																								
POCT FOB CARD RESULT 1		Negative																																									
POCT FOB CARD RESULT 2		Negative																																									
POCT FOB CARD RESULT 3		Negative																																									

- Lot # can be manually entered or scanned.
 - If value of Fail is selected in Internal QC, the system will default Abnormal in the Flags column as well as the Abnormal field under Results. **NOTE:** When correcting Results, you will need to manually clear the Abnormal flags, they will not clear automatically.
- Result Status defaults to final for all Enter/Edit Results. Select Accept on the bottom right.